


# PACE® Model 300

THREE CHAMBER TEMPORARY  
CARDIAC PACEMAKER

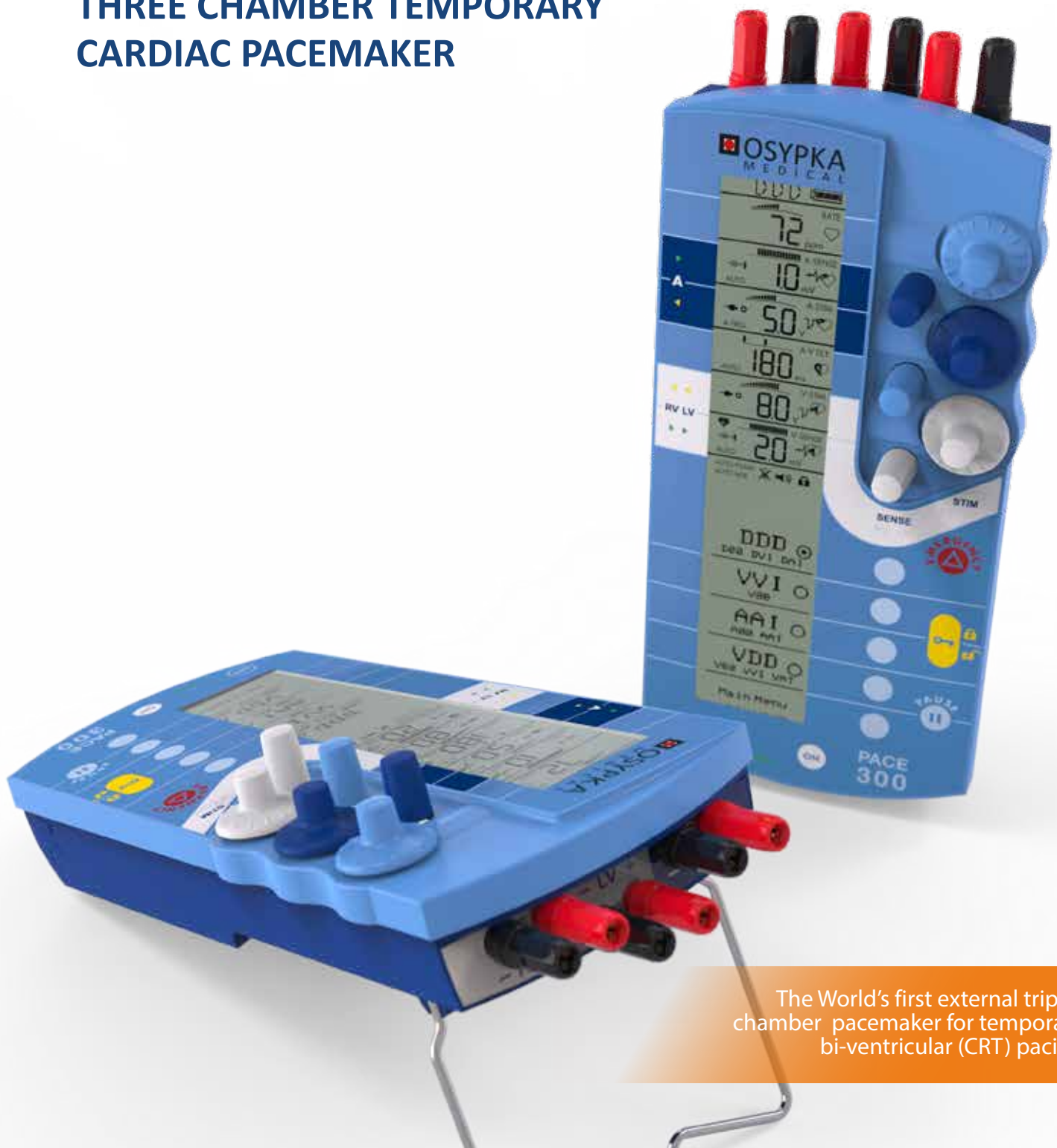

The World's first external triple-chamber pacemaker for temporary bi-ventricular (CRT) pacing

# Temporary Stimulation

Please refer to  
User's Manual for  
indications of use!

## PACE® Model 300 Highlights

- Three dedicated pacing terminals
- Resynchronisation in patients with low ejection fraction
- Constant voltage output for reliable cardiac capture (similar to permanent pacemakers)
- Variable setting of AV- and VV-Delay
- Auto Sense function automatically tracks P/R wave amplitudes and adjust atrial and ventricular sensitivities accordingly
- P and R Wave amplitudes measured upon pressing PAUSE key
- Quick and easy change of battery: no need to worry about battery orientation (one single 9 V battery)
- Lifetime service



P/N 2005

## Specifications

## PACE® Model 300

Modes	DDD (BV, R, L), VVI (BV, R, L), AAI, VDD (BV, R, L) D00 (BV, R, L), V00 (BV, R, L), A00, DVI (BV, R, L), DAI (BV, R, L), VAT (BV, R, L), AAT, DDD+AT, DAT (BV, R, L), DDI (BV, R, L)	
Basic Pacing Rates	30 – 220 ppm	
Rapid Atrial Pacing Rates	70 – 1,000 ppm	
Output Amplitude	0.1 – 18 V	
Pulse Width	0.05 – 1.5 ms	
V-V Delay	5 – 80 ms	
Sensitivity Threshold	Atrium : 0.2 – 20 mV	Ventricle: 1.0 – 20 mV
Refractory Period	250 ms	
Dimensions	21.7 cm x 9.6 cm x 5.1 cm (8.9" x 3.8" x 2.0")	
Weight incl. Battery	490 g (17.3 ounces)	
Battery Type	9V Type 6LP3146 Alkaline	
Battery Life	Alkaline: Up to 5 days (72 ppm, 8.0 V)	

Products and options may not be available in every region or country. Please contact [mail@osypkamed.com](mailto:mail@osypkamed.com) for availability in your region and country and for further information.

### Osypka Medical GmbH

Albert-Einstein-Strasse 3  
12489 Berlin, Germany  
Phone: +49 (30) 6392 8300  
Fax: +49 (30) 6392 8301  
E-Mail: [mail@osypkamed.com](mailto:mail@osypkamed.com)  
[www.osypkamed.com](http://www.osypkamed.com)

### United States of America: Osypka Medical, Inc.

dba. Cardiotronic  
La Jolla, CA 92037, USA  
Phone: +1 (858) 454 0600  
Fax: +1 (858) 454 0640  
E-Mail: [mail@cardiotronic.net](mailto:mail@cardiotronic.net)  
[www.cardiotronic.net](http://www.cardiotronic.net)

PACE and the logos are  
trademarks of Osypka Medical.